



**BIO MED  
WASH**®  
Bio-Logic Aqua®




FOR SKIN & EYES

**FIRST AID**

---

**STERILE WASH**

---

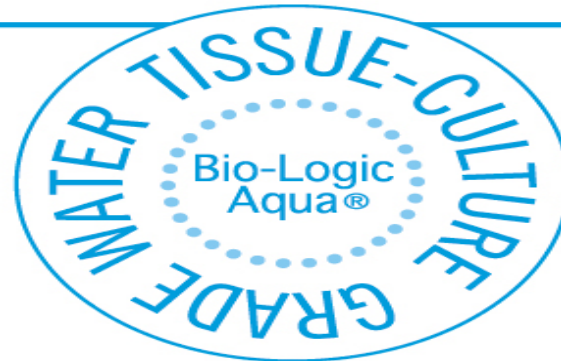



## Bio Med Wash is Unique

- **Personal Emergency Body or Eye Wash**
- Sterile Tissue Culture Grade Water w/ pH 6.4
- –ideal mineral nutrient content
- Saline-free, propellant-free & non-allergenic
- No impurities or contaminants quickly adsorbed by
- damaged tissue
- Trigger activated - Delivers fine spray which allows focused application in any direction or upside down!
- Replaces old & out-dated bottled saline solution
- Complies with American National Standard Z358.1-1998 as a Personal Eyewash





## **Bio Med Wash is inexpensive**

- Pennies per day cost of ownership
- Over 120 sprays per 7-oz can
- 4 year shelf life Vs 2 yr. for Saline & Buffered solutions
- Easy to deploy & use - Grab - snap - spray 2-4 seconds!
- Once opened does not need to be disposed of – stays sterile to the end!
- If frozen can be thawed and will not damage can
- Stations ready to mount - come with mounting plans or adhesive- Grab-N-Go bag for vehicles!



## **Bio Med Wash can be demonstrated**

- Great for safety meetings
  - Very few employees have ever seen an actual eye wash demo!
  - True attention grabber!
- Wonderful as corporate safety award!
  - Great for home, car, camp, boat and kids sporting events
- Personal spray for health concerns
  - Great for Computer dry eye & fatigue
  - Ready when and where you are!
- Very portable - assortment of sizes available



## Bio Med Wash size & types

- 3 convenient size applications (7-oz, 3-oz & 1-oz)
- 2 wall mount systems
  - 2-ea. 7-oz units with wall mount emergency station  
Self adhering or screw mount
  - 15-ea. 7-oz units with gravity dispenser system  
Designed to allow individual cans to be dispensed as needed
- 12 or 15 unit grab-and-go soft side bag  
Developed for utilities and remote sites
- Truck or wall mount for 7-oz  
Ideal for company vehicles, fork lifts, I-beams, labs & work cubicles
- Personal Nylon carry pouch 3-oz



## For Use By or In

- **Can be deployed anywhere you want or need additional worker protection**
  - Use in **conjunction** with plumbed & deployed eye wash units
- Ideal for Laboratories, vehicles, personal carry protection, tool boxes, first aid stations & remote work areas!
- Emergency Medical Technicians
- Plant Safety Programs
- PPE kits for Emergency Responders
- Police Officers and Correctional Officers

FEATURE	BENEFIT
<b>Available in three (3) Sizes</b>	1oz unit fits in first aid kits, 3 oz unit has a personal carry pouch, and 7oz units have truck mount, wall mount and 12 & 15 unit dispenser or carry pack.
<b>Sterile Tissue Grade Water</b>	BIO MED WASH® is light, saline free tissue-culture grade water proven to absorb more deeply and beneficially than other waters, including saline solution. Saline solution contains 140 times more dissolved solute than BIO MED WASH®. Although dissolved solutes in minute quantities are beneficial for healing, the more solutes water contains, the less the water is able to dissolve and wash impurities (saline solution is recommended for washing deep intramuscular cuts). Non-Saline.
<b>6.4 pH</b>	Slightly acidic, protects against eye irritation.
<b>Four (4) year shelf life</b>	Most personal eye wash units have only two year shelf life.
<b>Bag-In-Can Technology</b>	No Fluorocarbons, no gas release into atmosphere. Nitrogen bladder. Cutting edge "bag-in-can" technology improves sanitation, avoids chemical propellants and preservatives, and increases the product's useful life by insulating the water against extreme heat and cold.
<b>Trigger Activator</b>	Allows control usage of product. Multiple uses as product always remains sterile. The BIO MED WASH® trigger actuator emits Bio-Logic Aqua® tissue-culture grade water as a spray, making it unnecessary to touch the wound.
<b>Two (2) second activation and extended use</b>	Immediate ability to use. No covers to break open. 7oz unit contains over 100 single sprays or over 60 seconds of continuous use. Typical saline solutions are wrapped with cellophane shrink wrap. Once opens contents can drain in less than 10 seconds.
<b>Hydrotherapy and First Aid</b>	Bio Med WASH may be utilized for <b>eye irritation, burn management, wound management, tear gas &amp; pepper spray and to help rinse chemical exposure</b> . BIO MED WASH® all-natural, sterile first aid wash for cuts and lacerations, is the first product that replaces saline solution for washing abrasions and lacerations (saline is still preferred for deep intramuscular cuts). The product is invaluable in the initial washing of skin that has come in contact with irritating chemicals such as gasoline, pepper spray, tear gas or plant toxins; and the washing of eyes irritated by foreign particles, heavy smoke, dust or perspiration.

<b>FEATURE</b>	<b>BENEFIT</b>
<b>Contact lens use</b>	Will not wash out contact lens, help moisturize eyes.
<b>Protecting Your Eyes in an Emergency</b>	Quick and simple eye care information for fire-fighters, law enforcement, military and rescue personnel that could save vision - or lives!
<b>Clear Vision Can Save a Life</b>	It is a fact: Military personnel, law enforcement, fire-fighters, homeland security and emergency rescue workers are not taught enough about protecting their own vision during an emergency. The need is enormous because these well-trained and brave individuals are often faced with situations where "Clear Vision Can Save a Life™" – Not only their own but the lives of those around them.
<b>Vision Threats</b>	During a disaster or emergency, threats to vision may come from multiple sources: Smoke, chemical and explosive fumes, chemical weapons, tear gas, pepper spray, perspiration, dust, wind, heat, solar radiation, cold, stress, fatigue, foreign particles, etc. Vision threats may also originate from the more mundane: Building and window insulation, cleaning fluids, wet paint, auto exhaust, slicing onions. A plastic laundry basket left too close to the fireplace could create an emergency just as dangerous to life and vision as an oilfield explosion. When threat becomes reality, the resultant vision challenges can range from eyestrain to burning or itching eyes, to excessive tearing, to blurred vision. These challenges may be mild (sub-acute), temporarily debilitating or they could result in permanent visual impairment or loss of life.
<b>Fumes, smoke, tear gas, pepper spray, airborne chemicals, Chemical or thermal burns to eyes, eyelids, or skin</b>	Airborne chemicals create discomfort by alerting the pH balance (acidity/alkalinity) and osmolarity (moisture attracting and retaining ability), of tear film. Should the irritants build up to where discomfort becomes extreme irrigate with BIO MED WASH® by spraying directly into eyes from a distance of 8-12", followed by gentle blotting with a handkerchief or tissue. Protective eyewear or headgear also helps protect eyes from airborne irritants.
<b>BIO MED WASH® &amp; OC exposure</b>	BIO MED WASH® will probably not remove all OC ("oleoresin capsicum" or pepper spray) from skin or eyes. The product provides a sterile wash to start the decontamination process, and can also help reduce the risk of secondary exposure for someone attending a sprayed individual. Note: OC should be washed off with great care and according to established first aid procedures. Contaminated areas should not be touched or scrubbed.



# Bio Med Wash Pictures



Wall/Truck Mount



Soft Sided Bag  
-15 Pack-



3oz Holster



Eye Wash Magazine  
- Wall Dispenser -





**BIO MED**  
**WASH**®  
Bio-Logic Aqua®

